



The establishment of in vitro culture and drug screening systems for a newly isolated strain of *Trypanosoma equiperdum*


著者 (英)	Suganuma Keisuke, Yamasaki Shino, Molefe Nthatisi Innocentia, Musinguzi Peter Simon, Davaasuren Batdorj, Mossaad Ehab, Narantsatsral Sandagdorj, Battur Banzragch, Battsetseg Badgar, Inoue Noboru
journal or publication title	International Journal for Parasitology: Drugs and Drug Resistance
volume	7
number	2
page range	200-205
year	2017
URL	http://id.nii.ac.jp/1588/00001239/





Contents lists available at ScienceDirect

International Journal for Parasitology: Drugs and Drug Resistance

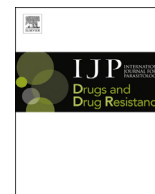journal homepage: www.elsevier.com/locate/ijpddr

The establishment of *in vitro* culture and drug screening systems for a newly isolated strain of *Trypanosoma equiperdum*



Keisuke Suganuma ^{a, b, *}, Shino Yamasaki ^b, Nthatisi Innocentia Molefe ^b,
Peter Simon Musinguzi ^b, Batdorj Davaasuren ^{b, c}, Ehab Mossaad ^{b, d},
Sandagdorj Narantsatsral ^c, Banzragch Battur ^c, Badgar Battsetseg ^c, Noboru Inoue ^e

^a Research Center for Global Agromedicine, Obihiro University of Agriculture and Veterinary Medicine, Inada, Obihiro, Hokkaido, Japan

^b National Research Center for Protozoan Diseases, Obihiro University of Agriculture and Veterinary Medicine, Inada, Obihiro, Hokkaido, Japan

^c Institute of Veterinary Medicine, Laboratory of Molecular Genetics, Mongolian University of Life Sciences, Zaisan 17024, Ulaanbaatar, Mongolia

^d Department of Pathology, Parasitology and Microbiology, Collage of Veterinary Medicine, Sudan University of Science and Technology, P.O. Box 204, Khartoum, Sudan

^e Obihiro University of Agriculture and Veterinary Medicine, Inada, Obihiro, Hokkaido, Japan

ARTICLE INFO

Article history:

Received 11 January 2017

Accepted 11 April 2017

Available online 13 April 2017

Keywords:

Colorimetric assay

Drug screening system

Liquid culture

Luciferase assay

Trypanosoma equiperdum

ABSTRACT

Dourine is caused by *Trypanosoma equiperdum* via coitus with an infected horse. Although dourine is distributed in Equidae worldwide and is listed as an internationally important animal disease by the World Organization for Animal Health (OIE), no effective treatment strategies have been established. In addition, there are no reports on drug discovery, because no drug screening system exists for this parasite. A new *T. equiperdum* strain was recently isolated from the genital organ of a stallion that showed typical symptoms of dourine. In the present study, we adapted *T. equiperdum* IVM-t1 from soft agarose media to HMI-9 liquid media to develop a drug screening assay for *T. equiperdum*. An intracellular ATP-based luciferase assay using CellTiter-Glo reagent and an intracellular dehydrogenase activity-based colorimetric assay using WTS-8 tetrazolium salt (CCK-8 reagent) were used in order to examine the trypanocidal effects of each compound. In addition, the IC₅₀ values of 4 reference trypanocidal compounds (pentamidine, diminazene, suramin and melarsomine) were evaluated and compared using established assays. The IC₅₀ values of these reference compounds corresponded well to previous studies involving other strains of *T. equiperdum*. The luciferase assay would be suitable for the mass screening of chemical libraries against *T. equiperdum* because it allows for the simple and rapid-evaluation of the trypanocidal activities of test compounds, while a simple, inexpensive colorimetric assay will be applicable in developing countries for the evaluation of the drug sensitivity of epidemic trypanosome strains.

© 2017 The Authors. Published by Elsevier Ltd on behalf of Australian Society for Parasitology. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

1. Introduction

Dourine is caused by *Trypanosoma equiperdum* of the subgenus *Trypanozoon*. Unlike other trypanosomes, dourine is transmitted through coitus with an infected horse rather than insect vectors (Brun et al., 1998). Thus, dourine has been found to be distributed

worldwide, beyond the tsetse belt in sub-Saharan Africa (Brun et al., 1998). Dourine is listed as an internationally important animal disease by the World Organization for Animal Health (OIE). Evidence from *in vitro* drug sensitivity assays using laboratory strains of *T. equiperdum* indicates that suramin, diminazene, quinapyramine and Cymelarsan[®] (melarsomine) are effective against

* Corresponding author. Research Center for Global Agromedicine, National Research Center for Protozoan Diseases, OIE Reference Laboratory for Surra, Obihiro University of Agriculture and Veterinary Medicine, Inada, Obihiro, Hokkaido 080-8555, Japan. Tel.: +81 155 49 5642; fax: +81 155 49 5643.

E-mail addresses: k.suganuma@obihiro.ac.jp (K. Suganuma), yamasakishino0824@yahoo.co.jp (S. Yamasaki), nthatisimolefe@gmail.com (N.I. Molefe), pmusinguzi@yahoo.com (P.S. Musinguzi), davlag_mgl@yahoo.com (B. Davaasuren), ehabmssd7@gmail.com (E. Mossaad), naran69@gmail.com (S. Narantsatsral), bat912b@yahoo.com (B. Battur), bata07@gmail.com (B. Battsetseg), ircpmi@obihiro.ac.jp (N. Inoue).

T. equiperdum (Zhang et al., 1991; Brun and Lun, 1994; Gillingwater et al., 2007). Melarsomine was found to be effective against the *T. equiperdum* OVI strain in mice, and the *T. equiperdum* Dodola strain in mice and horses (Hagos et al., 2010; Habte et al., 2015). However, this has not been established as an effective treatment strategy for dourine-infected horses. Moreover, diminazene- and melarsomine-resistant *T. equiperdum* strains have already been reported (Gillingwater et al., 2007). Thus, an *in vitro* drug screening assay for *T. equiperdum* is strongly needed for efforts to seek novel trypanocidal compounds both for evaluating drug sensitivity and for establishing an effective treatment strategy for dourine.

The direct determination of the trypanosome count using a cell counting chamber is a simple method for evaluating trypanocidal activity. This assay is relatively cheap and simple to perform; however, it is time consuming to evaluate large numbers of compounds using this method (Fumarola et al., 2004). The incorporation of a radioactive nucleotide assay has been widely used to analyze the drug sensitivity of trypanosomes *in vitro* (Brun and Lun, 1994). After the incubation of the trypanosomes with the drugs and ^3H thymidine, the signals of the radioisotope in the parasite are counted using a liquid scintillation counter. This test has the disadvantage of involving radioactive nucleotides and cumbersome steps. A fluorescence-based screening system using Alamar blue™ has also been developed for an HTS assay for trypanosomes (Sykes and Avery, 2009a). This assay requires expensive equipment, such as a fluorometer for the 96-well plate and is therefore difficult to apply in developing countries, which are most affected by trypanosomiasis (Fumarola et al., 2004). Intracellular ATP-based luciferase assays have also been utilized as HTS assays against *T. congolense*, *T. brucei* and *Leishmania* spp. (Sykes and Avery, 2009b; De Muylder et al., 2011; Sukanuma et al., 2014). The kits, such as CellTiter-Glo®, are simple, rapid and efficient, with signal detection taking only few minutes. Intracellular ATP-based luciferase assays are more suitable for the rapid-evaluation of trypanocidal activity than Alamar blue™ based fluorescence assays, because the incubation time after the addition of the CellTiter-Glo® reagents is short. However, the expensive equipment and reagents that are necessary for this assay also constrain its application in developing countries. Colorimetric assays using tetrazolium salts, such as MTT, XTT and WST, have also been established for *T. brucei* and *Leishmania* spp. (Lu et al., 2013; Ginouves et al., 2014). These methods are based on a color reaction (through the reductive reaction of intracellular dehydrogenase activity), which leads to the formation of a colored chemical product, formazan. The results of the color reaction are measured using a relatively inexpensive and versatile microplate plate reader (ELISA plate reader). In addition, the reagents are affordable in comparison to fluorescence and luciferase assays (Fumarola et al., 2004). Thus, these colorimetric assays are more suitable than other assay systems for evaluating the drug sensitivity of trypanosomes in developing countries.

We recently isolated a new strain of *T. equiperdum*, *T. equiperdum* IVM-t1 from a horse in Mongolia that was confirmed to have dourine based on clinical symptoms, a positive PCR result, the serological diagnosis, and the microscopic observation of trypanosomes in a genital swab (Sukanuma et al., 2016). This trypanosome strain was well adapted and propagated on the soft agarose media. In the present study, we successfully adapted *T. equiperdum* IVM-t1 that had been cultured in soft agarose media to culture in HMI-9 liquid media in order to develop a drug screening assay for *T. equiperdum*. In addition, the IC₅₀ values of four trypanocidal drugs were evaluated and compared to the results obtained from established assays.

2. Materials and methods

2.1. *In vitro* cultivation

The *T. equiperdum* IVM-t1 strain, which was isolated in Mongolia in 2015, was used in the present study (Sukanuma et al., 2016). Soft agarose culture-adapted *T. equiperdum* IVM-t1 was maintained by continuous sub-culturing once a week at 37 °C in 5% CO₂ using soft agarose media (HMI-9 (Hirumi and Hirumi, 1991) with 0.8% low gelling agarose (wt/vol) [Type VII, Sigma-Aldrich Japan, Tokyo, Japan]).

2.2. Adaptation to liquid media

Soft agarose culture-adapted *T. equiperdum* IVM-t1 parasites were transferred into HMI-9 liquid media supplemented with 0, 0.01, 0.1, 1 and 10% FcEM (vol/vol) (Nissan Chemical Industries, Ltd., Tokyo, Japan) at a parasite concentration of 1×10^4 cells/mL. The number of trypanosomes in each cell culture flask was manually counted every 24 h using a cell counting chamber. After optimizing the FcEM concentration in HMI-9, the liquid culture-adapted *T. equiperdum* IVM-t1 strain was maintained by replacing the entire culture supernatant with fresh HMI-9 supplemented with 1% FcEM every other day. A one-way ANOVA and Tukey's multiple comparison test were used for the statistical analyses of the doubling time, the daily parasite concentrations and the maximum parasite concentrations with different concentrations of FcEM. The statistical analyses were performed using the GraphPad PRISM 5 software program (GraphPad Software, Inc., CA, USA). The doubling time was calculated by counting the cells in the logarithmic growth phase of growth and using the following equation: Doubling time = $(t_2 - t_1) \times \log(2) / \log(q_2/q_1)$. Two measurements were performed: the initial count (q_1) at time (t_1) and the resultant parasite concentration after 24 h of incubation (q_2 , t_2).

2.3. The optimization of the luciferase assay

To optimize the efficacy of the parasite detection by the luciferase assay using the CellTiter-Glo® Luminescent Cell Viability Assay reagent (Promega Japan, Tokyo, Japan), trypanosomes were prepared by two serial dilutions from 2×10^6 – 1950 cells/mL, and 100 μL aliquots were dispensed into 96-well cell culture plates in duplicate. The luminescence was measured using 12.5, 25, 50 and 100 μL of CellTiter-Glo® reagents at various time intervals after 2 min mixing, in order to evaluate the linearity between the luminescence and the trypanosome concentration.

In order to optimize the initial trypanosome concentration in the 96-well cell culture plates, 100 μL of *T. equiperdum* at a density of 5×10^4 , 2.5×10^4 , 1×10^4 and 5×10^3 cells/mL was inoculated on a 96-well cell culture plate and the cells were manually counted every 24 h after inoculation.

Moreover, to optimize the concentration of DMSO (which was used as a solvent) in the culture medium, the inhibitory effect of DMSO was evaluated. One hundred microliters (1×10^4 cells/mL) of *T. equiperdum* were cultivated for 72 h in HMI-9 supplemented with 1% FcEM and 0, 0.125, 0.25, 0.5, 1, 1.25, 5 and 10% of DMSO. Subsequently, the survival rate of the parasite in each concentration was evaluated by a luciferase assay using 25 μL of CellTiter-Glo® reagent. Student's *t*-test was used to analyze the differences between the control group (0% DMSO) and DMSO-treated groups.

2.4. The assessment of the trypanocidal activity with optimized assay conditions

Fifty microliters of various concentrations of reference

compounds in HMI-9 supplemented with 1% FCEM and 0.5% DMSO were added to a Nunc™ MicroWell™ 96-well optical bottom plate (Thermo Fisher Scientific K.K., Yokohama, Japan) for the luciferase assay, or a 96-well cell culture plate (Corning™ Costar™ Flat Bottom Cell Culture Plates; Thermo Fisher Scientific) for the colorimetric and manual counting assays. Subsequently, 50 μL of liquid culture-adapted- or soft agarose culture-adapted *T. equiperdum* IVM-t1 in HMI-9 supplemented with 1% FCEM at a density of 2×10^4 cells/mL were aliquoted into each well and mixed by gentle pipetting in the 96-well cell culture plates. The final concentration of trypanosomes used in the assay was 1×10^4 cells/mL in 100 μL of HMI-9 supplemented with 1% FCEM and 0.25% DMSO per well. The trypanosomes were then incubated for 72 h. To evaluate the trypanocidal activity in the luciferase assay, 25 μL of CellTiter-Glo® was added and mixed for 2 min using an MS3 basic plate shaker (IKA® JAPAN K.K., Osaka, Japan). Subsequently, the plates were kept for 5 min in room temperature. The luminescence was then measured using a GloMax®-Multi + Detection System plate reader (Promega Japan). To evaluate the trypanocidal activity using the colorimetric assay, 10 μL of CCK-8 reagent was added, mixed for 2 min using an MS3 basic plate shaker and incubated in an incubator for 4 h in accordance with the manufacturer's protocol. Subsequently, the absorbencies in 450 nm were measured using an MTP-500 microplate reader (Corona Electric, Ibaraki, Japan). In addition, the trypanocidal activity was evaluated based on the number of living trypanosomes, which were manually counted using a cell counting chamber.

The trypanocidal activity of four commercial trypanocidal compounds (pentamidine [Sigma-Aldrich], diminazene aceturate [Sigma-Aldrich], suramin [Wako Pure Chemical Industries, Ltd., Osaka, Japan] and melarsomine [Merial, Lyon, France]) against *T. equiperdum* IVM-t1 was evaluated. These compounds were stored as 1 mg/mL (diminazene and suramin) or 40 $\mu\text{g}/\text{mL}$ (pentamidine and melarsomine) stock solutions in DMSO. Five-fold serial dilutions were made in duplicate from 2500 to 0.032 ng/mL (diminazene and suramin) and 100–0.00128 ng/mL (pentamidine and melarsomine) in HMI-9 with 0.25% DMSO. The IC_{50} values of each reference compound were calculated based on the results of the luciferase assay, the colorimetric assay and manual counting by plotting the % inhibition (0% inhibition = the luminescence, absorbance or trypanosome concentration of the trypanosome culture wells with 0.25% DMSO as a solvent control without reference compounds [negative control wells], 100% inhibition = the luminescence, absorbance or trypanosome concentration of trypanosome culture wells with 25 $\mu\text{g}/\text{mL}$ of diminazene [positive control wells]) in the GraphPad PRISM 5 software program using nonlinear regression (curve fit). A one-way ANOVA and Tukey's multiple comparison tests were used for the statistical analyses of the results. The statistical analyses were performed using the GraphPad PRISM 5 software program.

The reproducibility of these assays was validated using the coefficient of variation (%CV), signal/background ratio (S/B), signal/noise ratio (S/N) and Z'-factor values, which were calculated using the following equations.

$$Z': 1 - \frac{([3 \times \text{SD of negative control wells} + 3 \times \text{SD of positive control wells}])}{([\text{mean of negative control wells} - \text{mean of positive control wells}])}$$

as described by Sykes and Avery (2009a).

3. Results

3.1. The optimization of the *T. equiperdum* liquid culture conditions

Soft agarose media-adapted *T. equiperdum* IVM-t1 transferred into HMI-9 without FCEM reached a maximum cell concentration of $1.54 \pm 0.32 \times 10^6$ cells/mL (Fig. 1). The doubling time in the logarithmic growth phase in HMI-9 without FCEM was determined to be 16.36 ± 1.02 h. The doubling time was shortened by FCEM in a dose-dependent manner and the peak *T. equiperdum* concentration was increased by FCEM in a dose-dependent manner by supplementation with up to 1% FCEM. The doubling time of 1% FCEM culture (11.15 ± 0.23 h) was significantly shorter than that of other FCEM concentrations ($p < 0.05$). In addition, the peak trypanosome concentration of the 1% FCEM culture ($3.95 \pm 0.43 \times 10^6$ cells/mL; Day 4) was significantly higher than that of the other concentrations (Day 5, Fig. 1). These results suggested that HMI-9 with 1% FCEM was the most suitable media for the induction of liquid culture-adapted *T. equiperdum*.

3.2. The optimization of the in vitro trypanocidal activity evaluation system using CellTiter-Glo® reagent

The optimization of the volume of CellTiter-Glo® reagent and the time point for the measurement of luminescence was carried out by evaluating the correlation between the luminescence and the trypanosome concentration. This was achieved by adding various volumes of CellTiter-Glo® reagent and measuring the luminescence at different time points. The results showed that the coefficient of determination (R^2) values of 100 μL (the volume recommended by the manufacturer), 50 μL and 25 μL of CellTiter-Glo® reagent at 15 min after mixing were still ≥ 0.99 from the lowest concentration (1950 cells/mL) to the highest concentration (2×10^6 cells/mL) (Supplemental Fig. 1-A–C) Fig. 2-A shows the optimized condition, in which measurement was performed after the addition of 25 μL of CellTiter-Glo®, mixing for 2 min, and standing for 5 min at room temperature. From 7.5 min after mixing, the R^2 of 12.5 μL of the reagent was found to gradually decrease from 0.99 (Supplemental Figs. 1–D).

Next, the condition of the trypanosome culture in the 96-well plate that led to a final concentration of under 2×10^6 cells/mL in the logarithmic growth phase on day 3 was optimized by the cultivation of various concentrations of trypanosomes from day 0 to day 3 (Fig. 2-B). We cultured that were started with a concentration of $< 2.5 \times 10^4$ cells/mL yielded $< 2 \times 10^6$ cells/mL on day 3 while still in the logarithmic growth phase. The doubling time in the logarithmic phase and the peak of trypanosome concentration did not differ to a statistically significant extent among the initial trypanosome concentrations. Thus, a starting *T. equiperdum* concentration of 1×10^4 cells/mL in HMI-9 supplemented with 1%

%CV : Standard deviation (SD) of duplicated wells/Mean of duplicated wells.

S/B : Mean of negative control wells/Mean of positive control wells.

S/N : (Mean of negative control wells – Mean of positive control wells)/SD of positive control wells.

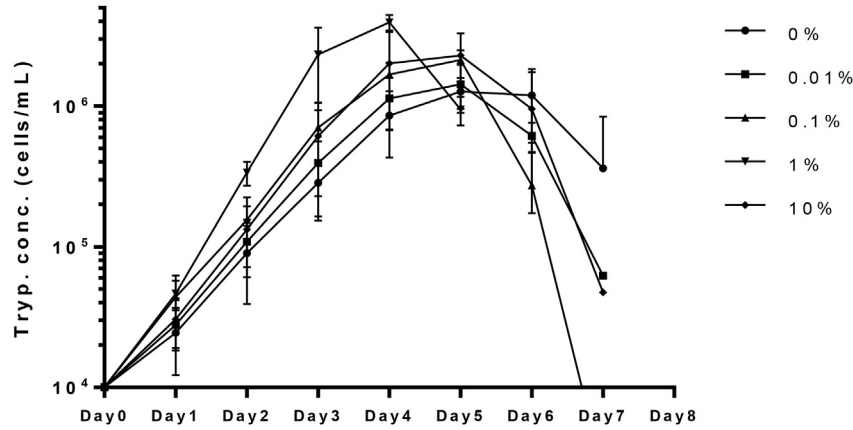


Fig. 1. Optimization of the FCEM concentration for liquid culture adaptation. Optimization of the FCEM concentration. Each plot shows the mean value from four independent experiments. Each column shows the mean value ± standard deviation.

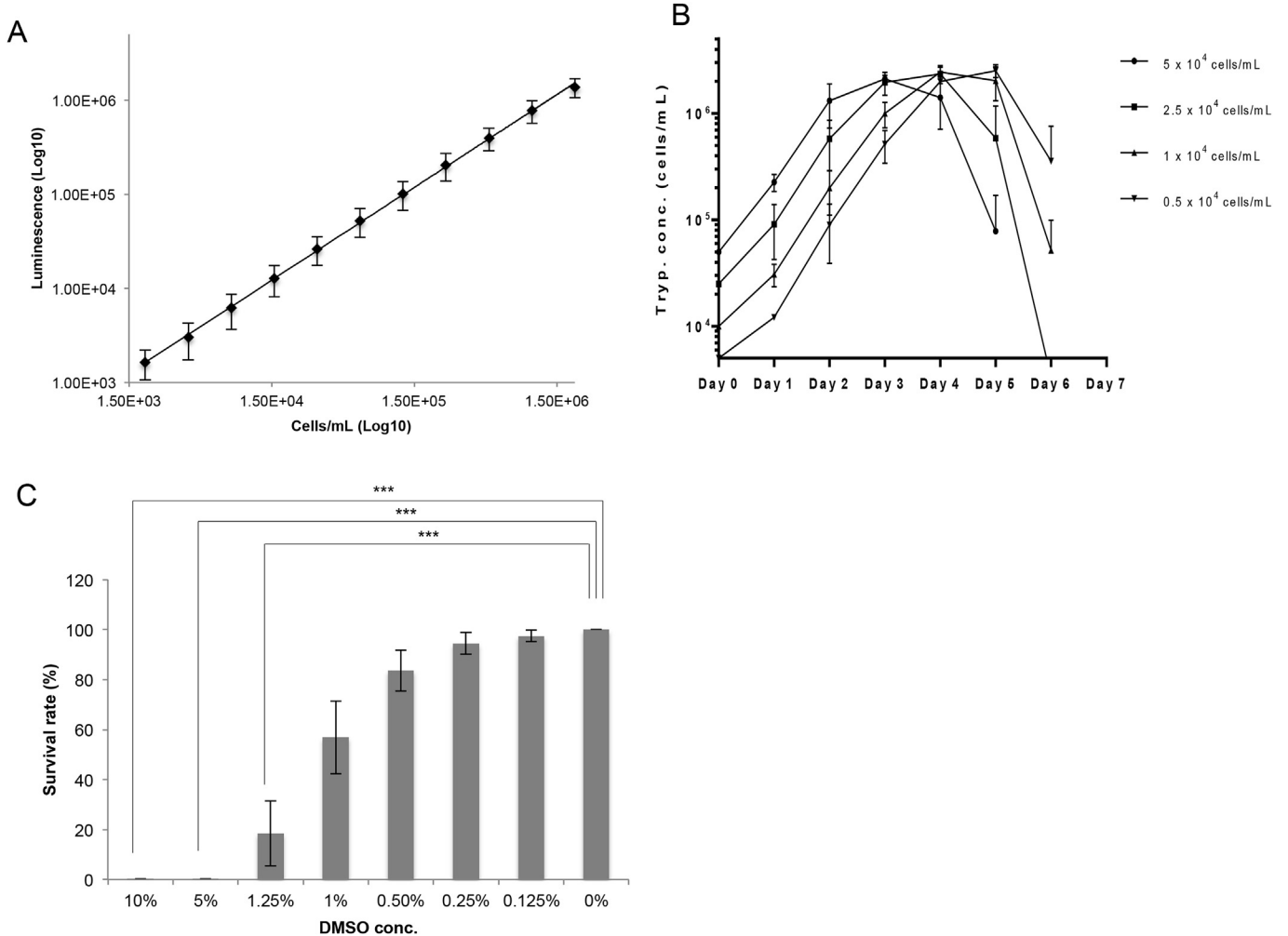


Fig. 2. The optimization of the conditions of the assay using CellTiter-Glo® reagent. (A) The correlation between the number of viable *T. equiperdum* IVM-t1 and luminescence was measured using the following optimized conditions: 25 µL of reagent at 5 min after mixing. All of the plots were calculated from four independent experiments and are shown as the mean luminescence signal ± standard deviation. The scale of the X and Y axes is log-10. (B) The optimization of the initial concentration for *T. equiperdum* IVM-t1 culture using 96-well cell culture plates. Each plot was calculated from three independent experiments and shows the mean value ± standard deviation. (C) The optimization of the DMSO concentration in the medium. All bars were calculated from three independent experiments and are shown as the mean survival rate ± standard deviation. Student's *t*-test was used to analyze the differences between 0% DMSO (negative control) and each of DMSO concentrations. ****p* < 0.001.

FCeM was determined to be the optimal condition for our trypanocidal screening assay using 96-well cell culture plates.

Finally, to optimize the DMSO (as a solvent of the reference trypanocidal compounds) concentration in the experiments, we evaluated the survival rate of *T. equiperdum* cultivated with various concentrations of DMSO. The supplementation of DMSO at concentrations of >1.0% significantly inhibited the growth of *T. equiperdum* ($p < 0.001$) (Fig. 2-C).

Based on these results, the economic efficiency of the reagents and the handling performance of assay, the optimal conditions for the determination of IC₅₀ values of the reference compounds were as follows: 25 µL of CellTiter-Glo[®] reagent, 0.25% DMSO, the measurement of luminescence after 2 min mixing and standing for another 5 min at room temperature.

3.3. The evaluation and comparison of the IC₅₀ values of the reference compounds against *T. equiperdum* IVM-t1

The IC₅₀ values against *T. equiperdum* were calculated by three different methods using *T. equiperdum* from two different origins: directly collected from soft agarose culture (soft agarose culture-adapted *T. equiperdum*); and collected from HMI-9 liquid culture supplemented with 1% FCeM (liquid culture-adapted *T. equiperdum*). The IC₅₀ value range of pentamidine was 0.65–0.86 ng/mL and no significant differences were observed among the methods and *T. equiperdum* from the two different origins (Table 1 and Supplemental Fig. 3). Furthermore, there were no significant differences in the IC₅₀ values of diminazene aceturate (IC₅₀ values, 4.16–6.16 ng/mL) (Table 1 and Supplemental Fig. 4). However, the IC₅₀ value of suramin, as evaluated by the colorimetric assays (26.75–35.49 ng/mL) was significantly lower than that determined by manual counting using liquid culture-adapted *T. equiperdum* (69.46 ± 5.19 ng/mL, $p < 0.05$) (Table 1 and Supplemental Fig. 5). Some of the IC₅₀ values of melarsomine also showed significant differences among the three methods and between the two different origins of *T. equiperdum* ($p < 0.05$) (Table 1 and Supplemental Fig. 6). The results of these drug screening assays were validated using the %CV, S/B, S/N and Z' values, as shown in Supplemental Table 1.

4. Discussion

Drug screening assay systems using liquid culture-adapted *T. equiperdum* are important for the development of novel trypanocidal compounds with *in vitro* activity against *T. equiperdum*. Bassarak et al. succeeded in establishing an efficient *in vitro* *T. equiperdum* culture system in a method that focused on the trypanosome metabolism by supplementation with a trace additive

(Bassarak et al., 2016). Recently, a new *T. equiperdum* strain, *T. equiperdum* IVM-t1, was established using a soft agarose culture system that was described in our previous study (Suganuma et al., 2016). This strain of *T. equiperdum* was characteristically well propagated with adherence on the surface of soft agarose, while the primary culture of this *T. equiperdum* strain could not be adapted to HMI-9 liquid media (Suganuma et al., 2016). In natural infection, *T. equiperdum* in the genital organs proliferates on the surface of the genital mucosa and invades tissues upon attachment (Brun et al., 1998). Thus, it was hypothesized that cellular adhesion is a key factor for adaptation and proliferation in *in vitro* culture. FCeM is a three-dimensional culture medium that contains FP001 polymer (Otsuji et al., 2014). Cells cultivated with media containing 0.015% FCeM (vol/vol) could maintain independence, with cells floating in a spheroid formation (Otsuji et al., 2014). Unlike other three-dimensional polymer, such as methylcellulose, the viscosity of the medium supplemented with FCeM is almost the same as water (Otsuji et al., 2014). Thus, we can handle trypanosome culture using HMI-9 supplemented with FCeM in the same way as normal HMI-9 medium for daily culture maintenance and further experiments. In FCeM supplemented medium, some *T. equiperdum* grow in a “cloudlike polymer-adherent colony” formation. The proliferation of *T. equiperdum* IVM-t1 with 1% FCeM supplementation was significantly better than that at lower concentrations of FCeM in HMI-9 media on days 3 and 4. In addition, the doubling time in the 1% FCeM-supplemented media was significantly shorter than in media supplemented with other concentrations of FCeM in the logarithmic growth phase. The peak trypanosome concentration with 1% FCeM supplementation (Day 4, $3.95 \pm 0.43 \times 10^6$ cells/mL) was significantly higher in comparison to the peak concentrations with the supplementation of other concentrations of FCeM on day 5 (0% FCeM, $1.54 \pm 0.32 \times 10^6$; 0.01% FCeM, $1.57 \pm 0.34 \times 10^6$; 0.1%, $2.19 \pm 0.23 \times 10^6$; and 10%, $2.23 \pm 0.93 \times 10^6$ cells/mL) (Fig. 1), and was almost the same as the cell concentration reported in a previous study (Bassarak et al., 2016). These results suggested that 1% FCeM supplementation in HMI-9 media was also the optimum media for the liquid culture of *T. equiperdum*.

The results of the optimization of conditions of the drug screening assay using CellTiter-Glo[®] reagents suggested that 1×10^4 cells/mL trypanosome as an initial culture concentration, and the measurement of luminescence using 25 µL of CellTiter-Glo[®] after 2 min mixing and 5 min standing at room temperature were suitable for the current study. Previous reports showed that the IC₅₀ values of diminazene, suramin and melarsomine (Cymelarsan[®]) against various strains of *T. equiperdum* (except for the drug-resistant OVI strain) were <0.1–19.2 ng/mL, 18.0–87.5 ng/mL and <0.1–2.0 ng/mL, respectively (Gillingwater et al., 2007; Brun and Lun, 1994). The IC₅₀ values of pentamidine and diminazene did

Table 1
The IC₅₀ values of the reference trypanocidal compounds.

		Pentamidine		Diminazene		Suramin		Melarsomine	
		Agarose	Liquid	Agarose	Liquid	Agarose	Liquid	Agarose	Liquid
Luciferase assay	ng/mL	0.79 ± 0.16	0.86 ± 0.11	6.16 ± 1.17	5.83 ± 1.47	43.96 ± 11.82	56.41 ± 17.37	^{a,c,d} 0.89 ± 0.01	^{c,e} 0.96 ± 0.10
	(nM)	(1.32 ± 0.27)	(1.46 ± 0.19)	(11.94 ± 2.27)	(11.30 ± 2.85)	(30.76 ± 8.27)	(39.47 ± 12.16)	(1.77 ± 0.22)	(1.91 ± 0.20)
Colorimetric assay	ng/mL	0.67 ± 0.13	0.72 ± 0.09	4.53 ± 0.36	5.04 ± 0.78	^a 26.75 ± 2.06	^{a,c} 35.49 ± 13.34	^b 0.49 ± 0.04	^{a,b} 0.71 ± 0.09
	(nM)	(1.13 ± 0.21)	(1.21 ± 0.17)	(8.78 ± 0.69)	(9.77 ± 1.52)	(18.72 ± 1.43)	(24.83 ± 9.34)	(0.98 ± 0.08)	(1.42 ± 0.18)
Manual Count	ng/mL	0.65 ± 0.10	0.76 ± 0.18	5.23 ± 1.21	4.16 ± 1.51	46.38 ± 17.43	^b 69.46 ± 5.19	^{a,b} 0.71 ± 0.15	^{a,c} 0.77 ± 0.12
	(nM)	(1.09 ± 0.17)	(1.29 ± 0.31)	(10.12 ± 2.34)	(8.06 ± 2.93)	(32.45 ± 12.20)	(48.60 ± 3.63)	(1.41 ± 0.29)	(1.52 ± 0.24)
Average	ng/mL	0.70 ± 0.15	0.76 ± 0.15	5.30 ± 1.25	5.01 ± 1.53	39.03 ± 15.69	53.79 ± 19.95	0.70 ± 0.20	0.81 ± 0.16
	(nM)	(1.18 ± 0.25)	(1.29 ± 0.25)	(10.29 ± 2.42)	(9.71 ± 2.98)	(27.31 ± 10.98)	(37.63 ± 13.96)	(1.39 ± 0.40)	(1.62 ± 0.31)

Each value shows the mean IC₅₀ value ± standard deviation from four independent experiments by ng/mL and nM.

Agarose, soft agarose culture adapted-*T. equiperdum* IVM-t1; Liquid, liquid culture-adapted *T. equiperdum* IVM-t1.

The different letters indicate that a significant difference was observed among the three methods and two origins in each reference drug ($p < 0.05$ Tukey's multiple comparison test).

not differ to a statistically significant extent among the methods and trypanosome origins (Table 1). Some statistically significant differences were observed among the IC₅₀ values of suramin and melarsomine that were obtained by different methods and in trypanosomes of the different origins; however, our IC₅₀ values mostly corresponded with those of previous reports (Table 1). In addition, the %CV, S/B, S/N and Z' values of the current assays were sufficient for a drug screening assay (Supplement Table 1). These results suggest that the established intracellular ATP-based luciferase assay using the CellTiter-Glo[®] reagent is suitable for high throughput screening (HTS) assays for *T. equiperdum*, if these optimized assay conditions can be applied in an automated screening system. Meanwhile, WTS-8 tetrazolium salt could be reduced and the color could react with the extracellular reductant in media. The HMI-9 media was supplemented with β-mercaptoethanol (as a reductant) to support trypanosome growth (Baltz et al., 1985). When following the manufacturer's protocol the linearity between the trypanosome concentration and absorbance was high ($R^2 \geq 0.99$) (Supplemental Fig. 2); however, a color reaction with the CCK-8 tetrazolium salt by β-mercaptoethanol caused high background noise and resulted in a low S/B value (1.61–1.90) (Supplemental Table 1). Moreover, if the test compounds have high reduction activity, the colorimetric assay using CCK-8 might not be applicable for an HTS assay because of the high background noise. Thus, it is strongly recommended that the trypanosomes be observed in 96-well cell culture plates by microscopy before evaluating the trypanocidal activity in an assay. Otherwise, the reference drug IC₅₀ values calculated by colorimetric assays mostly corresponded with those of previous reports (Table 1). Thus, in developing countries, the inexpensive colorimetric assay using CCK-8 might be more appropriate than the luciferase assay for evaluating the drug selectivity of epidemic trypanosomes.

In conclusion, we evaluated the methods of an HTS assay to determine the trypanocidal activity of different agents against *T. equiperdum*. Although there were some differences between the methods in terms of the IC₅₀ values of some of the reference compounds, both established assays were suitable for HTS assays against *T. equiperdum*. This offers a great advantage and an efficient tool for screening and determining novel dourine drugs for the treatment of dourine in future studies, not only in developed countries but also in developing countries.

Acknowledgements

We thank Ms. Yoko Matsushita for her excellent technical assistance. This study was financially supported by the Japan Society for the Promotion of Science (JSPS) KAKENHI Grant Number 16K18793 (Grants-in-Aid for Young Scientists [B]), and the AMED/JICA SATREPS. This research was partially supported by the “International Collaborative Research Program for Tackling the NTDs (Neglected Tropical Diseases) Challenges in African Countries” from

the Japan Agency for Medical Research and Development, AMED.

Appendix A. Supplementary data

Supplementary data related to this article can be found at <http://dx.doi.org/10.1016/j.ijppdr.2017.04.002>.

References

- Baltz, T., Baltz, D., Giroud, C., Crockett, J., 1985. Cultivation in a semi-defined medium of animal infective forms of *Trypanosoma brucei*, *T. equiperdum*, *T. evansi*, *T. rhodesiense* and *T. gambiense*. *EMBO J.* 4, 1273–1277.
- Bassarak, B., Moser, I., Menge, C., 2016. *In vitro* production of *Trypanosoma equiperdum* antigen and its evaluation for use in serodiagnosis of dourine. *Vet. Parasitol.* 223, 133–140.
- Brun, R., Hecker, H., Lun, Z.R., 1998. *Trypanosoma evansi* and *T. equiperdum*: distribution, biology, treatment and phylogenetic relationship (a review). *Vet. Parasitol.* 79, 95–107.
- Brun, R., Lun, Z.R., 1994. Drug sensitivity of Chinese *Trypanosoma evansi* and *Trypanosoma equiperdum* isolates. *Vet. Parasitol.* 52, 37–46.
- De Muylder, G., Ang, K.K., Chen, S., Arkin, M.R., Engel, J.C., McKerrow, J.H., 2011. A screen against *Leishmania* intracellular amastigotes: comparison to a promastigote screen and identification of a host cell-specific hit. *PLoS Negl. Trop. Dis.* 5, e1253.
- Fumarola, L., Spinelli, R., Brandonisio, O., 2004. *In vitro* assays for evaluation of drug activity against *Leishmania* spp. *Res. Microbiol.* 155, 224–230.
- Gillingwater, K., Buscher, P., Brun, R., 2007. Establishment of a panel of reference *Trypanosoma evansi* and *Trypanosoma equiperdum* strains for drug screening. *Vet. Parasitol.* 148, 114–121.
- Genouves, M., Carme, B., Couppie, P., Prevot, G., 2014. Comparison of tetrazolium salt assays for evaluation of drug activity against *Leishmania* spp. *J. Clin. Microbiol.* 52, 2131–2138.
- Habte, B., Bsrat, A., Ashenafi, H., Regassa, F., 2015. Efficacy of some trypanocidal drug against *Trypanosoma equiperdum* OVI in experimentally infected mice in debrezeit, Ethiopia. *EJBS* 7, 07–13.
- Hagos, A., Goddeeris, B., Yilkal, K., Alemu, T., Fikru, R., Yacob, H., Feseha, G., Claes, F., 2010. Efficacy of Cymelarsan[®] and Diminasan[®] against *Trypanosoma equiperdum* infections in mice and horses. *Vet. Parasitol.* 171, 200–206.
- Hirumi, H., Hirumi, K., 1991. *In vitro* cultivation of *Trypanosoma congolense* bloodstream forms in the absence of feeder cell layers. *Parasitology* 102 (Pt 2), 225–236.
- Lu, J., Vodnala, S.K., Gustavsson, A.L., Gustafsson, T.N., Sjöberg, B., Johansson, H.A., Kumar, S., Tjernberg, A., Engman, L., Rottenberg, M.E., Holmgren, A., 2013. Ebsulfur is a benzisothiazolone cytotoxic inhibitor targeting the trypanothione reductase of *Trypanosoma brucei*. *J. Biol. Chem.* 288, 27456–27468.
- Otsuji, T.G., Bin, J., Yoshimura, A., Tomura, M., Tateyama, D., Minami, I., Yoshikawa, Y., Aiba, K., Heuser, J.E., Nishino, T., Hasegawa, K., Nakatsuji, N., 2014. A 3D sphere culture system containing functional polymers for large-scale human pluripotent stem cell production. *Stem Cell Rep.* 2, 734–745.
- Suganuma, K., Allamanda, P., Hakimi, H., Zhou, M., Angeles, J.M., Kawazu, S.I., Inoue, N., 2014. Establishment of ATP-based luciferase viability assay in 96-well plate for *Trypanosoma congolense*. *J. Vet. Med. Sci.* 76, 1437–1441.
- Suganuma, K., Narantsatsral, S., Battur, B., Yamasaki, S., Otgonsuren, D., Musinguzi, S.P., Davaasuren, B., Battsetseg, B., Inoue, N., 2016. Isolation, cultivation and molecular characterization of a new *Trypanosoma equiperdum* strain in Mongolia. *Parasit. Vectors* 9, 481.
- Sykes, M.L., Avery, V.M., 2009a. Development of an Alamar Blue viability assay in 384-well format for high throughput whole cell screening of *Trypanosoma brucei brucei* bloodstream form strain 427. *Am. J. Trop. Med. Hyg.* 81, 665–674.
- Sykes, M.L., Avery, V.M., 2009b. A luciferase based viability assay for ATP detection in 384-well format for high throughput whole cell screening of *Trypanosoma brucei brucei* bloodstream form strain 427. *Parasit. Vectors* 2, 54.
- Zhang, Z.Q., Giroud, C., Baltz, T., 1991. *In vivo* and *in vitro* sensitivity of *Trypanosoma evansi* and *T. equiperdum* to diminazene, suramin, MelCy, quinapyramine and isometamidium. *Acta Trop.* 50, 101–110.